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
**Technical Data sheet for SURGISPON®**

|                                                    |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |                                                   |
|----------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------|
| <b>Company details</b>                             | <b>Name of Manufacturer:</b> Aegis Lifesciences Private Limited<br><b>Address:</b> 215 - 216, Mahagujarat Industrial Estate, At: Moraiya, PO: Changodar, TAL: Sanand, Ahmedabad – 382213 Gujarat, India.<br><b>Telephone number:</b> +91 9726903131<br><b>Fax:</b> +91 9727123131<br><b>E-mail Address:</b> admin@aegis-lifesciences.com                                                                                                                                                                                                                                                                                                                                                  |                                                   |
| <b>Products description and indication</b>         | <p>SURGISPON (Absorbable Gelatin Sponge) is a sterile, water insoluble gelatin absorbable sponge intended for haemostatic use by applying to a bleeding surface. It is biocompatible.</p> <p>SURGISPON is a surgical haemostatic sponge, manufactured from highly purified first grade gelatin material for use in various surgical procedures. When implanted in vivo and used in appropriate amounts, it is completely absorbed within &lt;4 weeks.</p> <p>SURGISPON can be effectively used in various surgeries for haemostasis, when control of capillary, venous and arteriolar bleeding by pressure, ligature and other conventional procedures is ineffective or impractical.</p> |                                                   |
| <b>Directive Class (under Directive 93/42/EEC)</b> | <b>III</b>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                | <b>Rule of classification: 17 (MDD 93/42/EEC)</b> |
| <b>Mode of action</b>                              | <p>SURGISPON is surface-mediated haemostatic devices. It arrests bleeding when applied directly to bleeding surfaces during surgery by providing a mechanical matrix that facilitates accumulation of blood and accelerating the clotting process of blood. The porous structure activates the thrombocytes at the moment blood comes in contact with the matrix of the sponge. This causes the thrombocytes to release a series of substances which promote their aggregation at the same time as their surfaces change character, thus enabling them to act as a catalyst for the formation of the fibrin.</p>                                                                          |                                                   |
| <b>Technical Specifications</b>                    | <p>Product is sterile, Off-white, odourless, biocompatible, Insoluble in water having water absorption capacity not less than 35 times of its weight of water and digestibility not more than 75 minutes, Residue on ignition is not more than 2.0%</p>                                                                                                                                                                                                                                                                                                                                                                                                                                   |                                                   |
| <b>Size/type:</b>                                  | <p>The Surgispon is manufactured in various sizes and shapes.</p>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |                                                   |
| <b>Composition / ingredients</b>                   | <p>Highly purified first extract grade porcine skin gelatin</p>                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |                                                   |

|                                                                       |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |
|-----------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Clinical safety</b>                                                | Product has history of safe use. Gelatin sponges are used since 1940 and till date they have undergone very little change. The US Pharmacopoeia contains monograph for Absorbable Gelatin Sponge. Many equivalent devices are available in market similar to SURGISPON. Clinical safety is further established by literature route of annex X of MDD directive.                                                                                                                  |
| <b>Biocompatibility</b>                                               | The device is biocompatible. Biocompatibility tests is carried out as per international standard EN ISO 10993-1:2009.                                                                                                                                                                                                                                                                                                                                                            |
| <b>Sterilization</b>                                                  | The product, SURGISPON is supplied in sterile form. Sterilization is carried out by gamma irradiation. It is determined that Gamma irradiation have no effect on the materials used in the product, its components and packaging. The minimum irradiation dose used for Surgispon sterilization is 25KGy and is carried out as per the method given in ISO 11137-2 and it does not adversely affect the quality, safety or performance of the product throughout its shelf life. |
| <b>Compatibility with drugs/biologics/liquids used during surgery</b> | The product is neutral and have no physio-chemical interaction with drugs/biologics generally used in surgical environment.                                                                                                                                                                                                                                                                                                                                                      |
| <b>Packaging description</b>                                          | SURGISPON is packed in a complete sterile barrier packaging system having high quality/medical grade quality packing material for all types of packing stages.                                                                                                                                                                                                                                                                                                                   |
| <b>Label</b>                                                          | As per ISO 15223-1:2016                                                                                                                                                                                                                                                                                                                                                                                                                                                          |
| <b>Instructions for use</b>                                           | See Attached IFU                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |
| <b>Precautions for handling or storage</b>                            | Handle with care, store at dry room at a temperature not more than 30°C.                                                                                                                                                                                                                                                                                                                                                                                                         |
| <b>Traceability</b>                                                   | GS1 Barcode                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |

Authorised Signatory

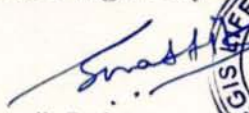Mr. Swastik Suthar  
Regulatory Affairs Assistant,  
Aegis Lifesciences Private Limited,  
Ahmedabad, India.  
15-10-2020